

Goodfellow Cambridge Limited


Ermine Business Park, Huntingdon PE29 6WR England Email: info@goodfellow.com

Telephone: +44 (0) 1480 424800

Fax: +44 (0) 1480 424900


This document will guide you through the easiest way to transfer graphene onto any compatible substrate

Graphene Transfer Procedure in 3 simple steps

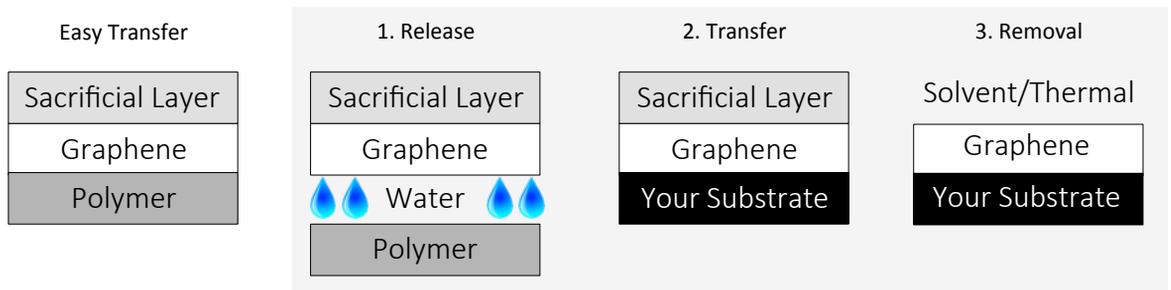

1. Release

Put the sample on deionized water slowly (no need to wet the sample first) while the sacrificial layer/graphene is detached from the support film. Once the sacrificial layer/graphene is floating remove the polymer film.

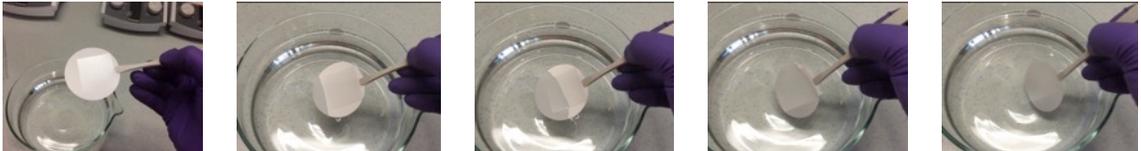

2. Transfer

The sacrificial layer/graphene will be floating in the water. To deposit it onto the desired substrate, introduce the substrate into the deionized water and fish the sacrificial layer/graphene from below. We recommend to tilt the substrate 45°. Take the sacrificial layer/graphene/substrate out and let it dry 30 minutes in air before store it under vacuum for at least 24h to avoid detachment of the graphene from your substrate.

3. Sacrificial layer Removal

To remove the sacrificial layer two methods can be used:

- **Solvents:** Dip the sacrificial layer/graphene/substrate in acetone for 1h, then into iso-propyl alcohol for another 1h and blow the sample with N₂ to dry it.
- **Thermal treatment:** Put the sacrificial layer/graphene/substrate into an oven and heat the sample at 450C in inert atmosphere for 2h.

Disclaimer: We believe that the information in this instruction is accurate and represents the best and most current information available to us. We make no representations or warranties either express or implied, regarding the suitability of the material for any purpose or the accuracy of the information contained within this document. Accordingly, we will not be responsible for damages resulting from use of or reliance upon this information.